


a new marker for histamine intolerance

HISTAMINE STOOL - URINE and WHOLE BLOOD Elisa Test

Enzyme Immunoassay for the quantitative determination of Histamine in stool, urine and whole blood. Histamine is quantitatively acylated. The competitive ELISA kit uses the microtiter plate format. The antigen is bound to the solid phase of the microtiter plate. Analyte concentrations of acylated standards, controls and samples and solid phase bound analyte concentrations compete for a fixed number of antiserum binding sites.

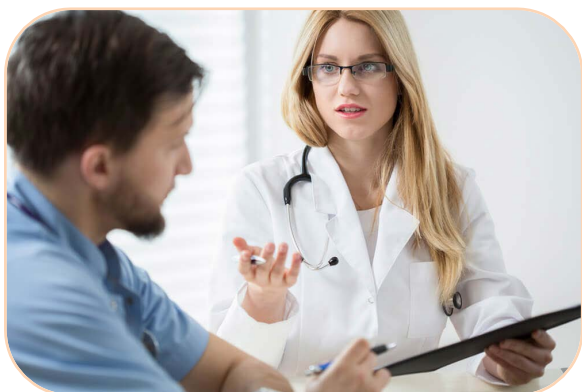

FEATURES

- Quantitative ELISA test
- Sample: Stool - Urine and Whole Blood
- ⌚ Incubation Time: 55 min
- Ideal for analysis of large number of samples
- Volume Sample: 25 µl of stabilized stool samples, urine or diluted whole blood
- Reagents, Standards and Controls : liquid ready to use
- Analytical Sensitivity: Stool 75 ng/g - Urine 0.9 ng/ml - Whole Blood 3.8 ng/ml
- 96 test/kit

nuovo biomarcatore per l'intolleranza all'istamina

ISTAMINA FECALE - URINE e SANGUE INTERO Elisa Test

Nuovo biomarcatore ELISA per la rilevazione quantitativa della concentrazione di Istamina su campioni fecali, urine e sangue intero. Il test consente di verificare in modo non invasivo, lo stato infiammatorio Intestinale e l'attivazione Immunitaria di reazione. Particolarmente indicato in caso di diarree ricorrenti, allergie diagnosticate o sospette, patologie autoimmuni sia a carico del tratto gastrointestinale che sistemiche.

SPECIFICHE

- Test quantitativo diretto su campioni fecali, urine e sangue intero.
- ⌚ Tempo del saggio: 55 min.
- Volume campione: 25 µl. di feci stabilizzate, urine e sangue intero.
- Reagenti, Calibratori e Controlli liquidi e pronti all'uso
- Sensibilità Analitica: Feci 75 ng/g - Urine 0.9 ng/ml - Sangue Intero 3.8 ng/ml
- 96 test/kit



EUROPE Ordering Information:

Histamine Stool - Urine - Whole Blood

Product Code: FDE-3000

Packaging: 96 test

Histamine Stool Collection tube

Product Code: FD-RO-40-GB

Packaging: 100 pcs.



Meridian Healthcare®

Meridian Healthcare srl

68 - Tremestieri Etneo

Catania, IT 95030

Phone: 095 -725 68 69 (0039)